
AGELOC® LUMISPA® ACNE SYSTEM 12-WEEK CLINICAL EFFICACY STUDY

Summarized from a third-party professional assessment performed in accordance with Good Clinical Practices and the Standard Operating Procedures of the Organization

©2019 NSE Products, Inc., Provo, UT

INTRODUCTION

Nu Skin created a novel skin care device that can improve skin appearance while cleansing the skin. It does so through Micropulse Oscillation technology—a pulsing, counter-rotation movement of 15 Hz. This results in a unique movement of the skin, which is an interplay between a silicone surface and specially created topicals. We assessed the benefits of this skin movement on the appearance of acne-prone skin throughout a 12-week clinical evaluation.


METHODS

Nu Skin Research and Development contracted with a third-party testing house in the United States to evaluate the efficacy of the ageLOC LumiSpa Acne System. Thirty male and female subjects with Fitzpatrick skin types I-IV and between 18–45 years of age with acne-prone skin were selected to use the ageLOC LumiSpa Gentle treatment head with the ageLOC LumiSpa Treatment Cleanser—Acne. Subjects were instructed to use the ageLOC LumiSpa System in place of their normal facial cleansing morning and night. Clinical grader and self-perception evaluations were conducted on an array of skin beauty characteristics at baseline and weeks one, two, four, eight, and 12.


RESULTS

Using the Investigator Global Acne Assessment (IGA), each subject's acne was evaluated by a clinical grader. The subjects came into the study with mild to moderate IGA assessments. Statistically significant improvements were noted in as early as one week. These improvements grew to a 71% over the course of the 12-week study. The average rating of the subjects' combined IGA Assessment at the end of the study was clear to almost clear compared to baseline scores. This system was an effective resolution of acne.


Investigator Global Acne Assessment

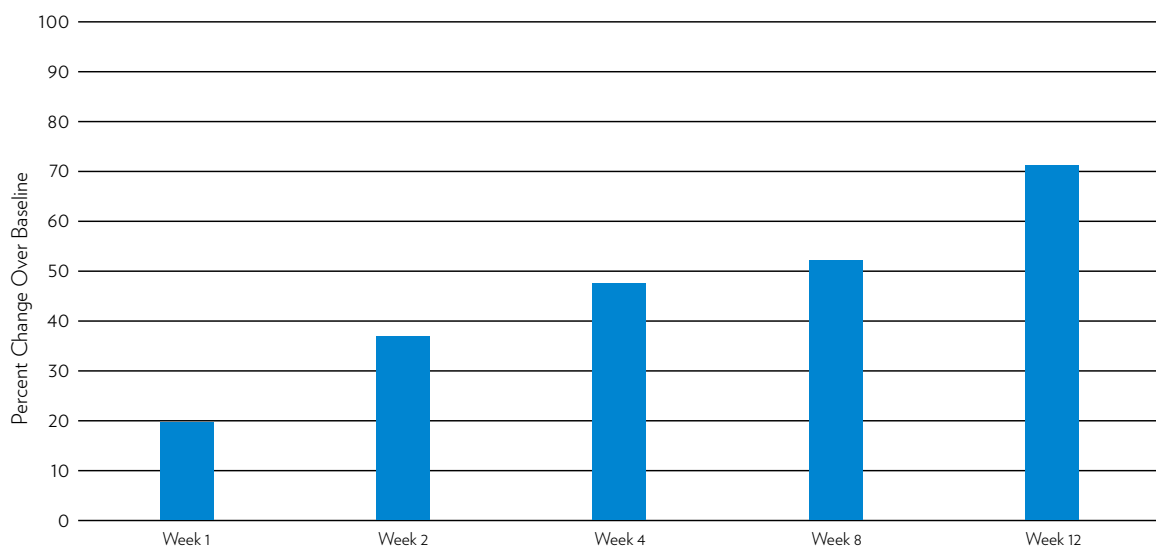

These same subjects were graded for hydration and cleansing ability. The clinical grader noted a 100% improvement in cleansing by week two which remained for the duration of the study. Similarly, these subjects were graded for a 100% improvement in hydration at week eight and twelve, illustrating the ability of the system to balance the skin.

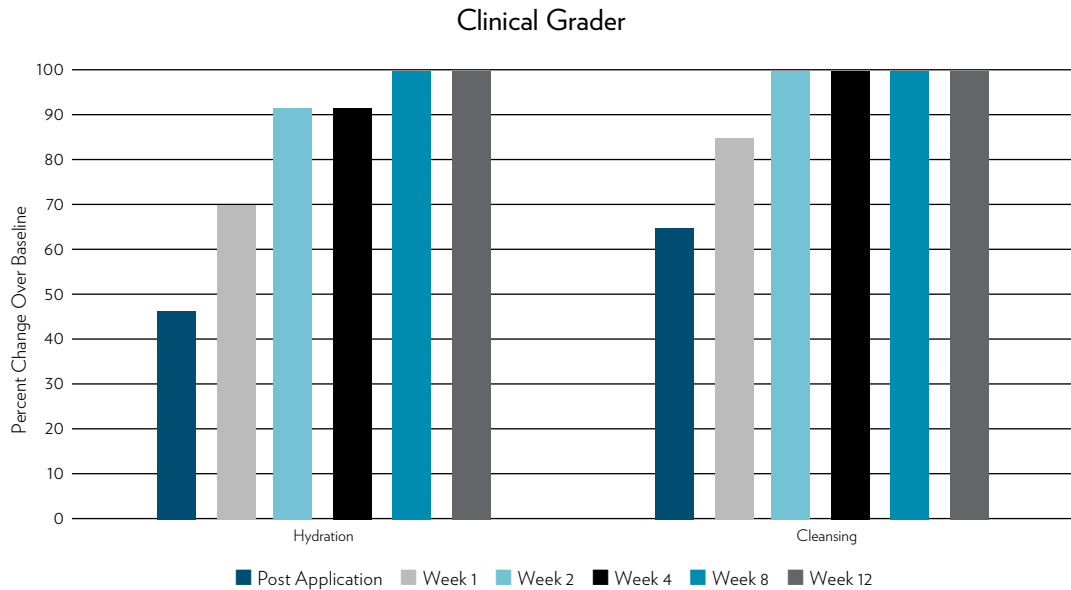

The clinical grader was also able to see visible improvements in skin smoothness, softness, pores, clarity and radiance.

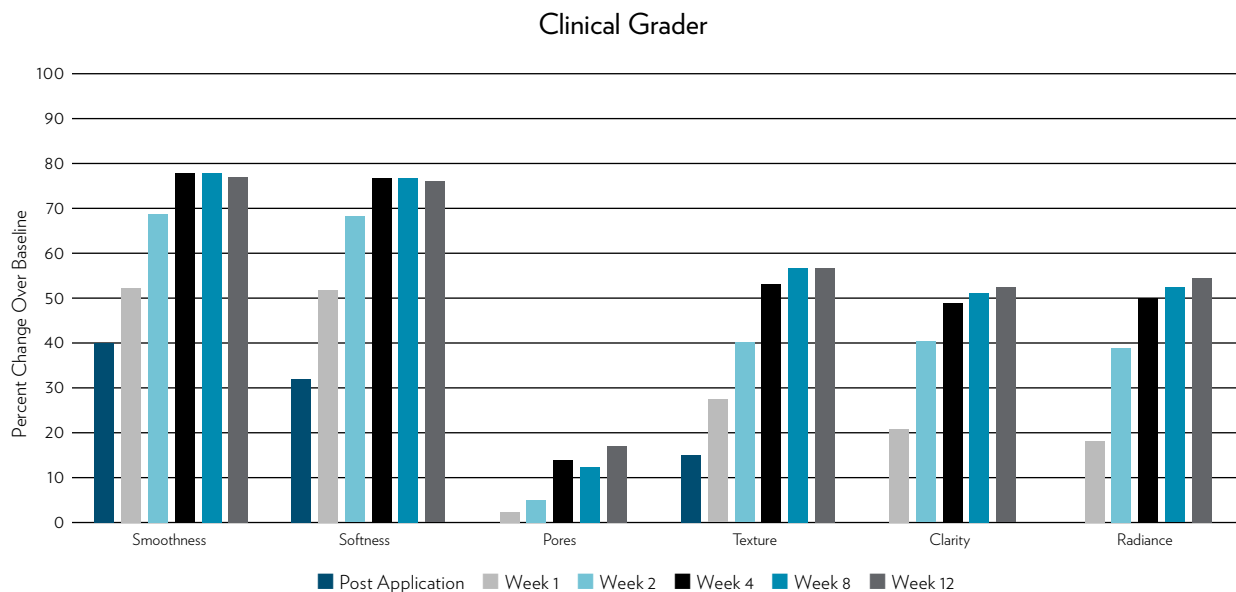

CONCLUSION

ageLOC® LumiSpa® Acne system, ageLOC LumiSpa with the Gentle head, and ageLOC LumiSpa Treatment Cleanser—Acne can help mitigate acne breakouts and help the skin look clearer.

